



DISPOSABLE CLEANING BRUSHES

FOR RELIABLE CLEANING

The disposable cleaning brushes from MICRO-TECH ensure a perfectly cleaned endoscope, thus reducing the risk of cross-contamination to a minimum. The brushes also protect your endoscope. Spherical beads on the brush ends prevent the endoscope's working channel

from being damaged during cleaning. If you wish to make your work even safer, more rapid and more efficient, the disposable double cleaning brushes with 2 heads are recommended. Two cleaning processes can therefore be performed simultaneously in a single working step.

SPECIFIC CHARACTERISTICS

- 2 brush heads for time-saving cleaning
- Protective beads on the brush heads
- Brushes made of resistant nylon
- Individually hygienically packed



Handle



Nylon brushes

CONICAL BRUSH HEADS

The product range has also been supplemented with disposable cleaning brushes with tapered brush heads. According to preference, brushes with constant and decreasing brush diameters are available. Conical brush heads are indicated by a T in the product number.




YOUR CHOICE: 1 OR 2 BRUSH HEADS

An ergonomic ring handle on the single head brush supports easy removal of the cleaning brush from the endoscope channel. Alternatively the use of a double headed cleaning brush can be recommended - a catheter with a brush at each end that is pulled through the channel thus allowing greater cleaning efficiency with its double brush action.



PRECISELY FITTING DIAMETER

Whether a bronchoscope or a therapeutic endoscope with a large working channel is involved – you can obtain ideal cleaning brushes for all endoscopes with diameters ranging from 3 mm to 11 mm.

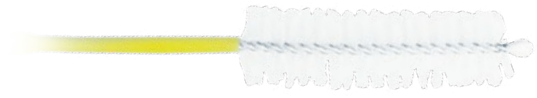

SPECIFICATIONS

| REF | Catheter Ø mm | Total length mm | Brushes Ø mm | Brush length mm |
|--|---------------|-----------------|-----------------|-----------------|
| CLEANING BRUSHES WITH ONE HEAD | | | | |
| CB18-0050/18 | 1.8 | 1800 | 5 | 20 |
| CB18-0050/23 | 1.8 | 2300 | 5 | 20 |
| CB18-T0050/23 | 1.8 | 2300 | 5-7 tapered | 20 |
| CLEANING BRUSHES WITH TWO HEADS | | | | |
| CB18-3030/18 | 1.8 | 1800 | 3/3 | 20/20 |
| CB18-5050/18 | 1.8 | 1800 | 5/6 | 20/20 |
| CB18-5050/23 | 1.8 | 2300 | 5/5 | 20/20 |
| CB18-T5050/23 | 1.8 | 2300 | 5-7/5-7 tapered | 20/20 |
| CB23-5050/18 | 2.3 | 1800 | 5/5 | 20/20 |
| CB23-5050/23 | 2.3 | 2300 | 5/5 | 20/20 |

Packaging unit: 100 pieces