

COMBINATION CLEANING BRUSH SOLUTIONS

THE 2 IN 1 SOLUTION FOR THE DAILY PRACTICE

The cleaning of both the endoscope channels and its air-water valve requires a different kind of brush. MICRO-TECH offers smart devices for this purpose that makes the daily

practice easier and more efficient. The dual cleaning brush and valve brush combination set made by MICRO-TECH is such a solution for the reprocessing of the endoscope.

SPECIFIC CHARACTERISTICS

- Double cleaning brush for air-water valve/control head
- Intensive and efficient cleaning
- Protective beads on the brush heads
- Brushes made of resistant nylon
- Individually hygienically packed sets available



Tapered



Protective bead

ONE BRUSH FITS ALL VALVES


MICRO-TECH offers the all inclusive solution for the thorough brushing of the air-water valve vents and ports. The double header brush combines a 40 mm long and 11 mm thick port brush with a 20 mm long and 5 mm thin vent brush. Both brushes are flexible and have protective distal beads to help prevent damage to the endoscope.



Double headed cleaning brush for air-water valve

COMBINATION DUAL FUNCTION BRUSH

A cleaning brush with a conical brush head for the endoscope channels and a large calibre and long brush head for the airwater valve- ports allow channel and valve-port brushing with a single cleaning brush. This saves time and reduces costs.

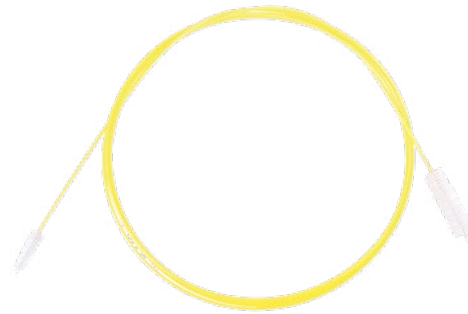

Combination cleaning brush

SMART AND EFFICIENT: CLEANING BRUSH SETS

The MICRO-TECH cleaning brush set includes all necessary brushes for a thorough and easy cleaning of the endoscope. It contains one cleaning brush for the endoscope channels and one double header cleaning brush for the valve vents

and ports. The sets are alternatively available with one single or one double headed cleaning endoscope channel brush. Each set is delivered individually hygienically packed.

SPECIFICATIONS

REF	Catheter Ø mm	Total length mm	Brushes Ø mm	Brush length mm
CLEANING BRUSH FOR VALVE VENT				
CB13-50120/01	1.3	160	5/12	20/40
COMBINATION CLEANING BRUSH				
CB18-T50110/23	1.8	2300	5-7 tapered/11	20/40
CLEANING BRUSH SET				
CB18-T0050/23-B				
CB18-T0050/23	1.8	2300	5-7 tapered	20
CB13-50120-01	1.3	-	5/12	20/40
CB18-T5050/23-B				
CB18-T5050/23	1.8	2300	5-7 tapered	20/20
CB13-50120/01	1.3	-	5/12	20/40

Packaging unit: 100 pieces